



Transfer with dignity
Quality Strength Versatile




PATIENT LIFTS



Porter 500

The Porter-500 is a mobile Transfer Lift designed and developed in collaboration with Health Care Professionals. The sleek, strong built lift, included with a complete range of functional slings, offers caregivers flexibility and versatility, for a great variety of Patient/Resident lifting with care.

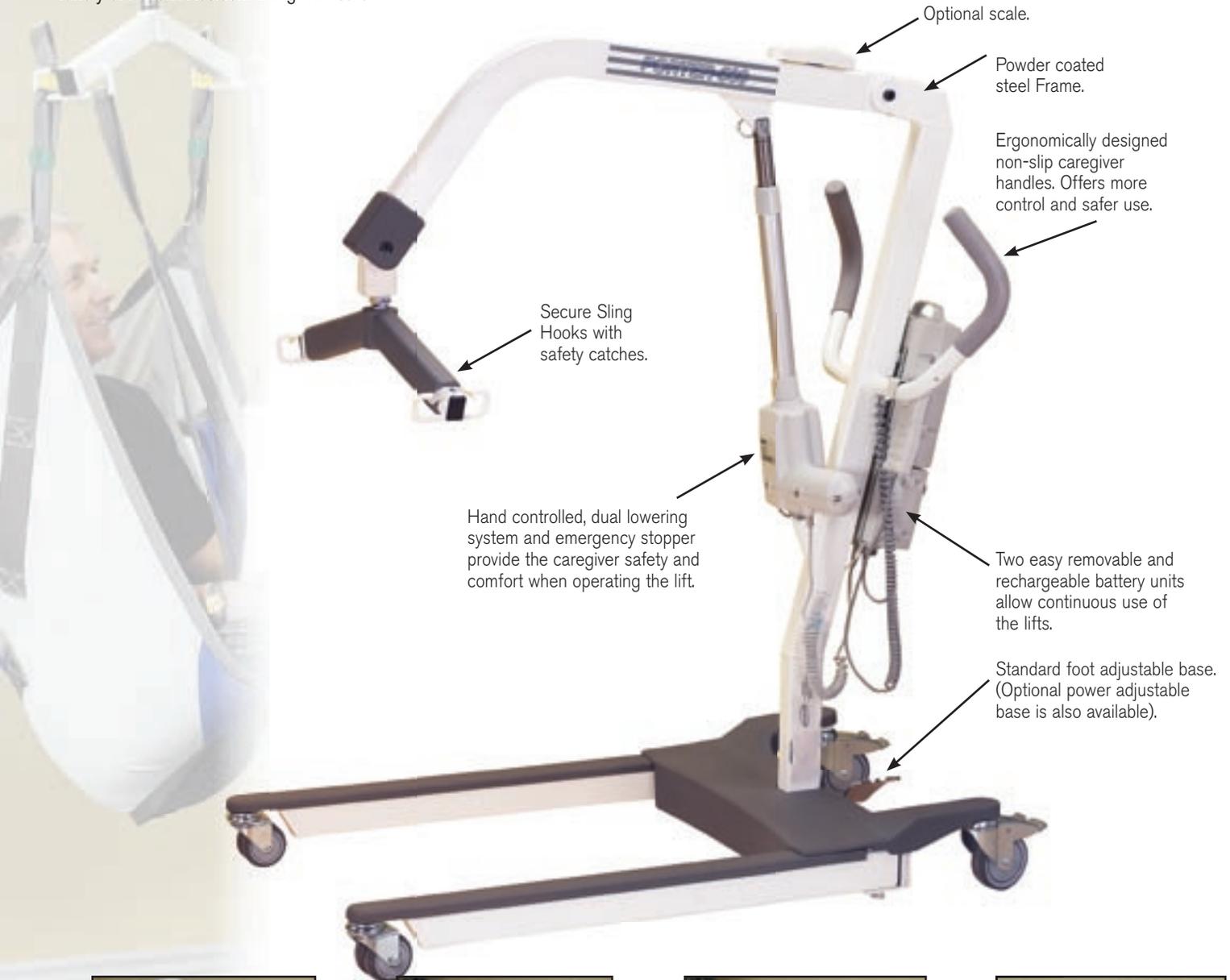

The Porter-500 is ideal for most lifting situations. With reduced patient/resident stress, the caregiver is able to provide secure, safe dignified lifts between beds, chairs, wheelchairs and toilets.



The Porter-500 is designed to provide minimal effort for manoeuvring in tight spaces while the spreading legs accommodate wide chairs and toilets.




The Porter-500 centres around the patient. The extra long boom provides more lifting height creating a more vertical lifting motion, and helps to prevent caregiver back injuries especially when lifting from the floor. Patients/Residents are more comfortable during transfers



Wall mounted Charging Station for batteries and adapter 120V, (230V adapter for Euro Plug is available).

Helper 400

The Helper-400, a versatile sit/stand lift, was designed using state of the art technology, developed in collaboration with health care professionals. The mobile, compact design offers a combination of lifting options, maximizing patient/resident safety and comfort.

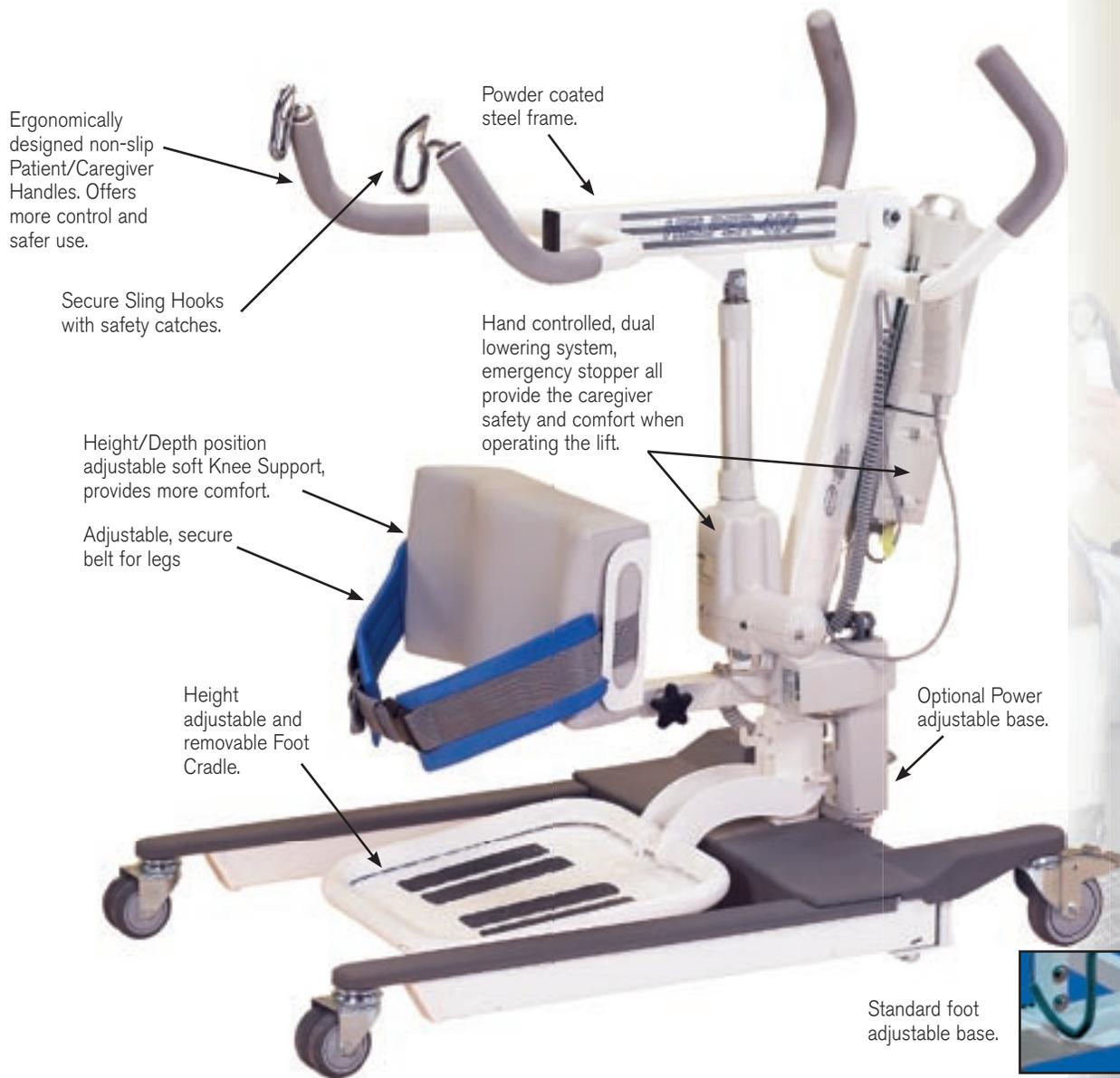

The Helper-400 offers caregivers quick, easy and safe patient/resident transferring between beds, chairs, wheelchairs, and toilets.



The Helper-400 reduces patient/resident and caregiver stress, allows more dignified dressing and toileting, and promotes secure, safe transferring. Compact design provides easy manoeuvring in tight spaces and spreading legs accommodate wide chairs and toilets.



Two easy removable and rechargeable battery units allow continuous use of the lift.




Wall mounted Charging Station for batteries and adaptor 120V. (230V adaptor for Euro Plug is available).

Wispa™ Lite Portable Hoist

The versatile moving and handling solution for the home, hospital and care environment.

The Wispa Lite is a portable hoist with a 200kg lifting capacity. It is the ideal moving and handling solution with value added functionality and a compact, lightweight design.

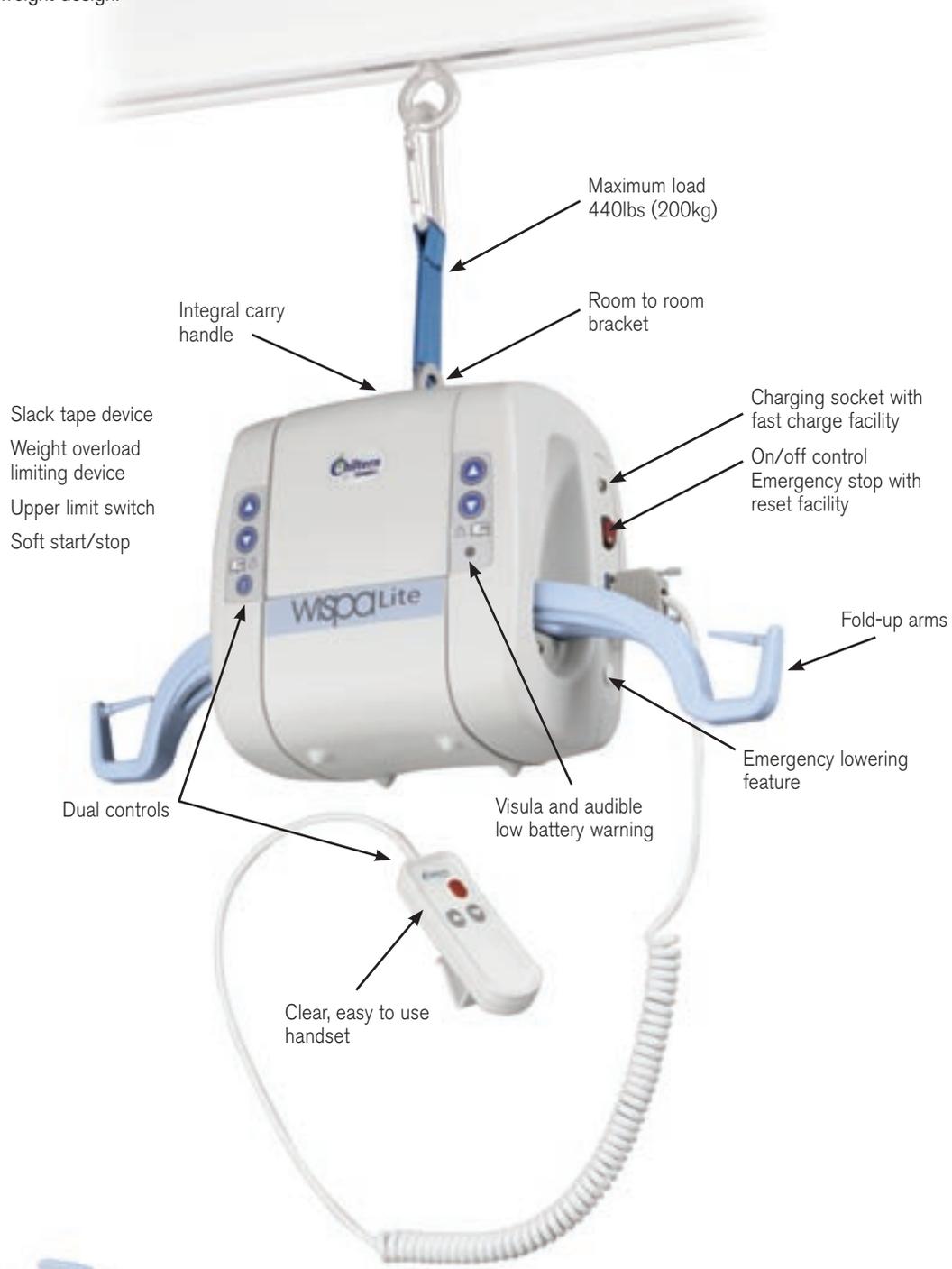

Accessories:

- Transfer Trolley
- Attachment handle
- Lanyard (length options)
- Comprehensive range of Slings, Standing Jackets and Accessories

Adaptable:

Easy to Transport

Lightweight construction combined with an integral carry handle enables the Wispa Lite to be easily moved between different locations, attached/removed from the ceiling track and repositioned wherever it is required.



Power Lift

Powerful, controlled lifting - Maximum load 440lbs (200kg).

Room to Room

The Wispa Lite can be moved effortlessly from room to room, without the need to reposition the client.

Easy to Connect to Track

With varying lengths of lanyard and attachment handle, caregivers can comfortably attach the Wispa Lite to ceiling track, whatever its height.

Varied Usage

Compatible with ceiling track for freestanding gantry systems offering controlled transfer from the bed, wheelchair or floor.

Caregiver & User Friendly:

Lightweight

The Wispa Lite is easy to carry, attach and remove from the ceiling track.



Handset

Clearly identifiable buttons, ergonomically designed.

Smooth Lifting

Soft start/stop lifting for added user comfort, safety and control.

Comfortable Lifting

Enhanced sling attachment design offers wider sling positioning to improve user comfort, whilst providing maximum head clearance.

Upper Limit

Switch

Prevents the lifting tape from overwinding into the case and up to the ceiling track.

Slack Tape Device

This mechanism prevents the lifting tape from winding out of the cassette when an obstruction is encountered during lowering, preventing jamming of the tape.

Emergency Lowering

In the event of hoist failure, Wispa Lite is fitted with an emergency lowering switch providing safe, controlled descent.

Dual Controls

Wispa Lite can be operated via controls on the cover or from the handset, thereby reducing risk of not lifting/lowering clients if the handset is lost or damaged.

Emergency stop with reset facility

Safe operation, should an emergency arise, providing reassurance for both user and caregiver. The hoist is easily reset via the handset or controls on the hoist.



Weight overload limiting device

If too much weight is applied, the Wispa Lite will only allow lowering.

Easy recharge facility

Visual and audible warnings indicate when the hoist should be recharged. With its integrated recharging feature, the Wispa Lite can simply be plugged into the nearest outlet.

Fast Charge

Wispa Lite can be fully charged within 2 hours.

Robust Design

Vigorously tested to ensure long product life expectancy.



Self-diagnostics

Audible (from hoist and handset) and visual signals (on hoist cover) indicate faults and battery status.



Flexible Track Options & Accessories


Turntable: Adding flexibility to track layout. Manual or powered model options. Points: Changes track direction. X-Y System: Efficient room coverage. Bends: Available in 15°, 30°, 45°, 90° sections.

Comprehensive Sling Range

Lightweight, flame retardant and durable material. Colour coded straps for quick and simple size recognition. Paediatric sizes and material available.

Bariatric Solutions - Wispa Twin Cassette

Safe working load: 998lbs (445kg). Free standing gantry system.

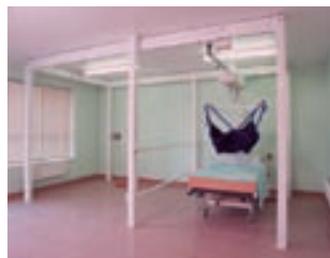

Fixed Ceiling Track Hoist - Wispa™ 100/200/300 Series

Maximum load: 100 series 286lbs (130kg)

200 series 440lbs (200kg)

300 series 598lbs (272kg)

Manual and powered travers options.

Optional Infrared control.

Service and Installation


Our nationwide team of trained engineers can quickly install your equipment, and service and maintain it to comply with LOLER regulations, on a 365 days/24 hour per day basis.

Warranty/guarantee

Wispa Lite comes with a two year guarantee (conditions apply).

Extended Warranty Programmes

Chiltern Invadex Super Care Warranty covers the full cost of breakdown (both parts and labor) for a period of up to five years. Please ask for details.

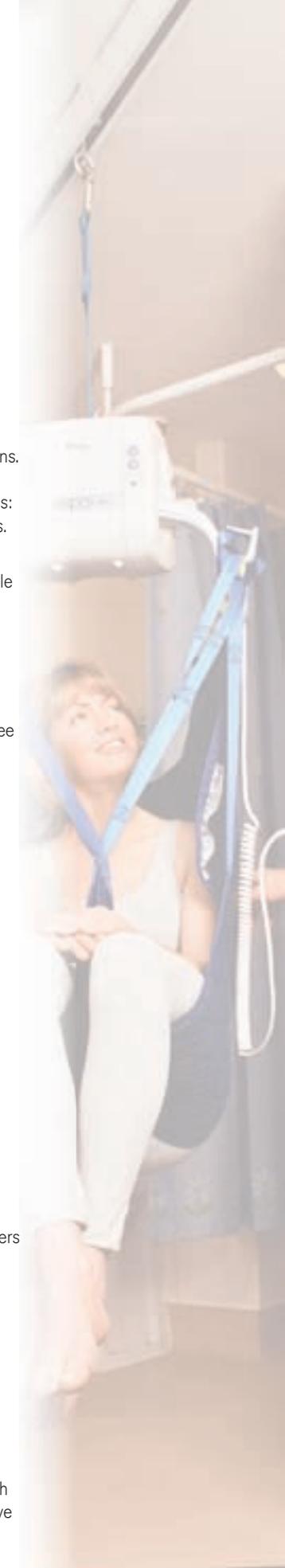

Wispa™ 200 and 300 Series Ceiling Mounted Hoist Range

Lifting solutions for the home, hospital and care environment.

One unique range, four different model options. Each intelligently designed to suit both user and caregiver. The Wispa Range is easy to operate with soft start and stop, the range allows safe, controlled transfer from a bed, chair, toilet or bathing area with precise fingertip control. Choose from manual traverse, powered traverse, or X-Y.



Powerful but Controlled: The 200 range lifts up to 440lbs (200kgs) or for the heavier patient the professional 300 model lifts up to 598lbs (272kgs). For added safety a limiting device is built-in to restrict the lift capability.



Finger Tip Control

Provided by the new waterproof handset (submersible up to 1 metre). Clearly marked for safe, efficient handling it features an assist button, which sounds an audible alarm in case of emergency.

Soft start/stop

No abrupt halt at the end of travel as the Wispa has a soft start and stop facility ensuring safe, controlled movement.

Variable Speed

The speed of the hoist travel can be adjusted, to suit the individual. For accurate positioning the variable speed facility is activated by the handset, allowing the hoist to glide effortlessly to the required position.

Economical but Simple to Recharge:

Always ready to go, Wispa is easy to charge. Simply return the hoist to the track end to start the charging process. An audible alarm sounds when the hoist needs recharging.

Dynamic Yet Safe: Responsive performance combines with a host of safety features, including:

Limiting Device

To stop the maximum load being exceeded. If the operator is distracted a limit switch prevents the patient from being raised too high or lowered too far.

Emergency Stop with Reset Facility

Should an emergency arise, this simple feature provides reassurance for both user and caregiver.

Emergency Lowering


In case of battery failure a lowering device is fitted.

Lift Motor has a Self Lock Facility


To prevent the user from being lowered unintentionally.



X-Y model option allows for complete room coverage.

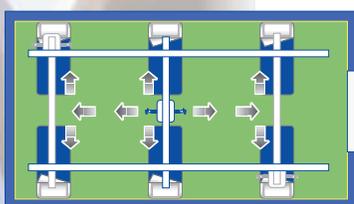


The "XY" system provides complete room coverage. For example, a user could be transferred to any bed in a hospital ward.

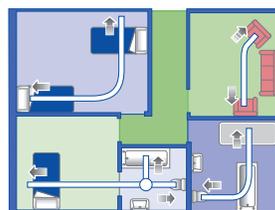

Moving between rooms is made simple. The Wispa allows easy transfer, as ceiling track can be fitted in most rooms throughout the house.

Wispa™ LiteTrack Gantry System

The Wispa LiteTrack gantry combines with the Lite Portable hoist to create the ideal moving and handling solution for the home, hospital and care environment. Totally freestanding, LiteTrack is ideal as a semi-permanent, temporary or emergency solution or where a fixed ceiling track installation or mobile hoist is inappropriate.

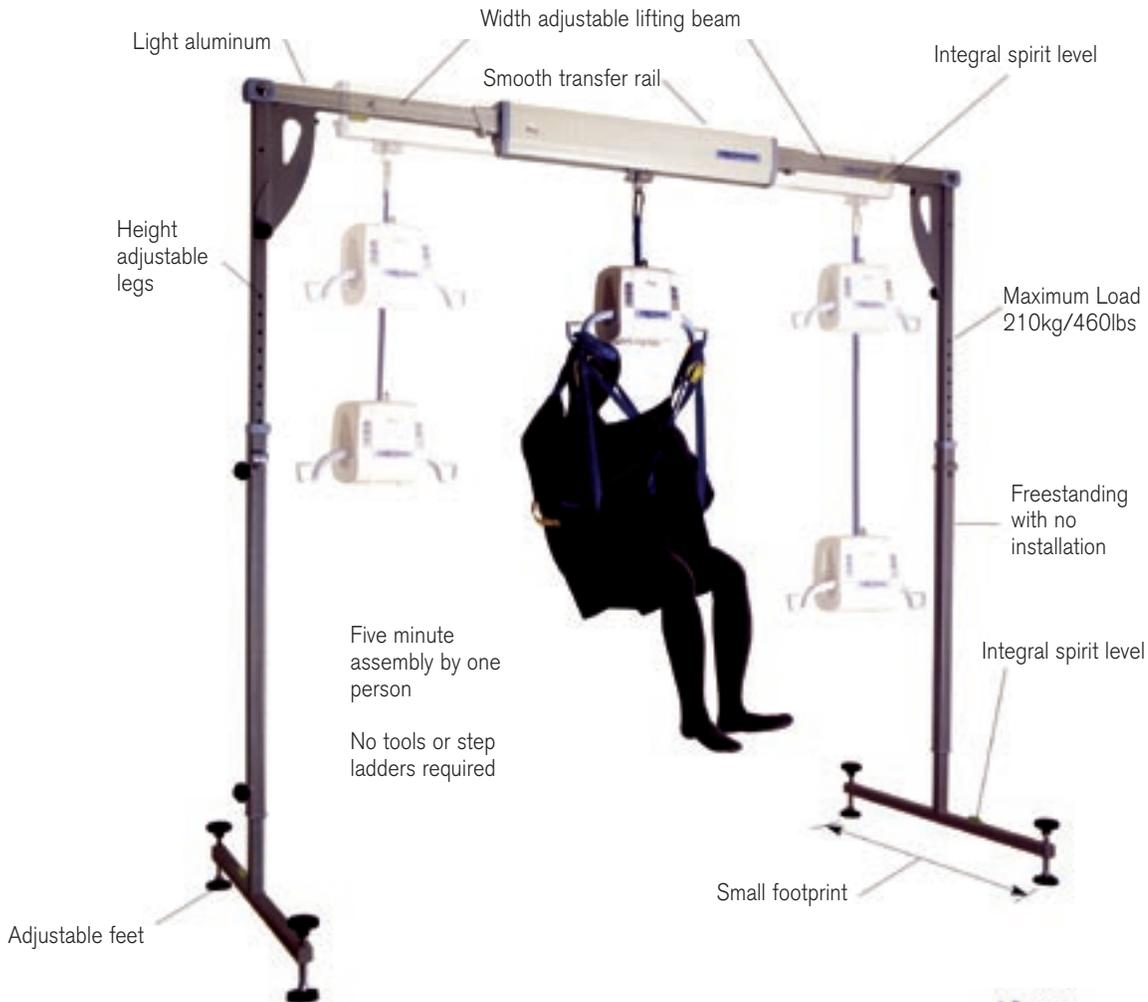

Suits the Majority of Clients

The LiteTrack gantry system complements the Wispa Lite portable hoist and patient slings with its maximum load of 460lbs (210kg).

Freestanding System/No Installation

Totally self supporting and independent of the building structure. Provides an alternative where the ceiling construction limits fixed track installation, floor construction, positioning issues or heavy/difficult clients limit the use of a mobile hoist. No installation requirement, electrical connection or fixings and can be removed in minutes. Operates with a battery powered hoist.

Easy Assembly

Transforms easily from a lightweight transportable system into a fully assembled gantry, in line with the requirements of the manual handling regulations. Unlike some gantry systems LiteTrack can be assembled or disassembled by one person within five minutes and without the use of tools or step ladders.

Lightweight Component Weights

The gantry is constructed from lightweight aluminum (the heaviest component part weighs only 24lb (11kg) with steel legs for added strength and an unobtrusive slim profile.

Compact and Transportable

For use in the community, the gantry has a small footprint to suit most room layouts and fits neatly into the trunk or back seat of a car, or for storage within emergency service vehicles.

Height and Width Adjustable

Ten pre-set height options at 50mm intervals to maximise the available ceiling height and three pre-set width options each side at 13-mm intervals to suit the lifting task and room layout.

Adjustable to Suit the Floor Surface Uneven floor surfaces are easily overcome with the adjustment wheels on the feet and an integral spirit level.



Compact and transportable



Ease of transfer



Wispa Lite portable hoist options